MitoFit Preprint Arch (2019) MiP/MitoEAGLE School Coimbra2019 2019-06-15 doi:10.26124/mitofit:ea19.MiPSchool.0005

Partitioning between cytochrome c oxidase and alternative oxidase studied by oxygen kinetics of dark respiration in Chlamydomonas reinhardtii: a microalgae model organism

<u>Marco Di Marcello</u>¹, Iglesias-Gonzalez J¹, Meszaros A^{1,2}, Haider M³, Gnaiger E^{1,2}, Huete-Ortega M¹

http://www.mitofit.org/index.php/Di Marcello 2019 MitoFit Preprint Arch



© 2019 Di Marcello *et al.* This is an Open Access extended abstract (not peerreviewed) distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original authors and source are credited. © remains with the authors, who have granted MitoFit an Open Access preprint license in perpetuity.

Editor MitoFit Preprint Archives: Gnaiger E

Introduction

Bioenergetics is the study of how living organisms acquire and transform energy to perform biological work. Energetic coupling between chloroplasts and mitochondria has been described in algae, demonstrating that a good functionality and interaction between both organelles is necessary to maintain metabolic integrity. High-resolution respirometry (HRR) is widely used to assess mitochondrial respiration and other bioenergetics parameters in the biomedical field of mitochondrial research and its clinical applications [1]. In our interdisciplinary study, we adapted the multimodal approach of the Oroboros O2k high-resolution respirometer to investigate algal bioenergetics for biotechnological purposes [2].

In contrast to mammalian cells, algal mitochondria possess alternative oxidases (AOX), which bypass electron transfer from the Q-junction through Complexes CIII and CIV [3]. Therefore, in algae we can distinguish between AOX-dependent and cytochrome c oxidase-dependent respiration through the Q-AOX and CIII-CIV pathways.

Material and methods

The microalgal model organism *Chlamydomonas reinhardtii* wild-type strain wt12 was grown at RT in Tris-Acetate-Phosphate (TAP) medium in a 16:8 h light:dark cycle. Oxygen flux, J_{02} , was monitored in wt12 living cells in the exponential growth phase at 25 °C in Oroboros O2k high-resolution respirometers excluding any light in the chambers. Substrate-uncoupler-inhibitor titration (SUIT) protocols were specifically developed to characterise activities of the Q-AOX pathway (SUIT-022 [4]) and CIII-CIV pathway (SUIT-023 O2 [5]). To quantify the contribution of the Q-AOX pathway to algal dark respiration, we studied the oxygen kinetics of (1) ROUTINE-respiration in TAP medium, (2) Q-AOX dependent

respiration after inhibition of CIV with 1 mM potassium cyanide (KCN), and (3) dependent respiration after inhibition of AOX with salicylhydroxamic acid (SHAM). Oxygen kinetics was obtained from aerobicanaerobic transitions with high time resolution at a data sampling interval of 0.2 s. p_{50} is the O₂ partial pressure, p_{O2} , at 50% of maximal respiration, J_{max} [6]. The p_{50} was calculated from hyperbolic fits using the Oroboros O2Kinetics software for automatic O₂ calibration, correction for zero O₂ signal drift, instrumental background O₂ flux and exponential time constant of the polarographic oxygen sensor [7]. A single shifted hyperbolic fit was used to fit J_{02} as a function of p_{02} in each aerobic-anaerobic transition.

Results and conclusions

 p_{50} ranged from 0.06 to 0.08 kPa for ROUTINE-respiration with an excellent fit by a first-order hyperbolic function. This oxygen affinity is comparable to that in small mammalian cells [8]. Upon inhibition of CIV with KCN, J_{02} was significantly impaired (Fig. 1A) and p_{50} increased three-fold up to 0.35 kPa (Fig. 2). No decline of J_{02} and p_{50} was observed relative to ROUTINE-respiration after inhibition of AOX with SHAM (Fig. 1B). In all cases, excellent fits of respiration as a function of oxygen pressure were obtained by a first-order hyperbolic function.

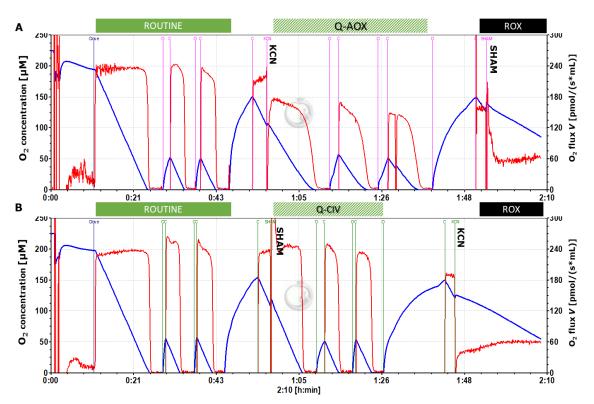


Figure 1. High-resolution respirometry for the study of dark respiration and O_2 kinetics with C. reinhardtii wt12. Representative O2k traces showing O_2 concentration and O_2 flux per chamber volume with repeated aerobic-anoxic transitions (O_2 kinetics) and re-oxygenations. A: Protocol SUIT-022: AOX-ce CN+SHAM. B: Protocol SUIT-023: AOX-ce SHAM+CN. Note the high technical reproducibility of ROUTINE-respiration in both protocols, and the identical and relatively high residual oxygen consumption, Rox, after titration of both inhibitors in both protocols.

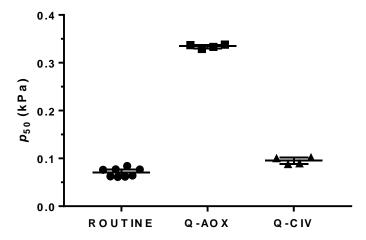


Figure 2. p_{50} values obtained for different metabolic pathways in *C. reinhardtii*. O₂ kinetic experiments were run in presence of the cytochrome c oxidase inhibitor potassium cyanide (AOX group) or the alternative oxidase inhibitor salicylhydroxamic acid (CIV group). The data represents n=8 technical replicates, N=2, median \pm IR.

If the potential contribution of the Q-AOX pathway in the ROUTINE-state would be compensated for by increased CIII-CIV pathway flux after addition of SHAM, then the mixed Q-AOX and CIII-CIV pathways would give rise to biphasic hyperbolic oxygen kinetics, with a contribution of the high-affinity CIII-CIV pathway and the low-affinity Q-AOX pathway. Taken together, our results provide evidence against a contribution of AOX to ROUTINE-dark respiration in wt12 cells under the presently applied culture conditions. Oxygen kinetics provides a fast and simple method for detection of AOX and CIV activities in dark respiration of living microalgal cells, extending our current understanding of the different O_2 affinities of the two pathways in these organisms and their possible effects on the bioenergetics and metabolism of the cells.



Supported by project NextGen-O2k, which has received funding from the European Union's Horizon 2020 research and innovation programme under grant agreement No 859770

References

- 1. Doerrier C, Garcia-Souza LF, Krumschnabel G, Wohlfarter Y, Mészáros AT, Gnaiger E (2018) High-Resolution FluoRespirometry and OXPHOS protocols for human cells, permeabilized fibers from small biopsies of muscle, and isolated mitochondria. Methods Mol Biol 1782:31-70. www.bioblast.at/index.php/Doerrier 2018 Methods Mol Biol
- 2. Pulz O & Gross W (2004) Valuable products from biotechnology of microalgae. App. Microbiol. Biotech. 65(6): 635-648.
- 3. Young L, Shiba T, Harada S, Kita K, Albury MS, Moore AL (2013) The alternative oxidases: simple oxidoreductase proteins with complex functions. Biochem Soc Trans 41:1305-11.
- 4. www.bioblast.at/index.php/SUIT-022 O2 ce D051
- 5. www.bioblast.at/index.php/SUIT-023 O2 ce D053
- 6. Gnaiger E (2001) Bioenergetics at low oxygen: dependence of respiration and phosphorylation on oxygen and adenosine diphosphate supply. Respir Physiol 128:277-97. www.bioblast.at/index.php/Gnaiger 2001 Respir Physiol

- 7. Meszaros AT, Haider M, Di Marcello M, Gnaiger E (2018) High-resolution mitochondrial oxygen kinetics as diagnostic tool in Complex IV impairments. Abstract Mitochondrial Medicine 2018 Hinxton UK.
 - www.bioblast.at/index.php/Meszaros 2018 Mt Med Hinxton
- 8. Scandurra FM, Gnaiger E (2010) Cell respiration under hypoxia: facts and artefacts in mitochondrial oxygen kinetics. Adv Exp Med Biol 662:7-25. www.bioblast.at/index.php/Scandurra 2010 Adv Exp Med Biol