

## Partitioning between cytochrome c oxidase and alternative oxidase studied by oxygen kinetics of dark respiration in *Chlamydomonas reinhardtii*: a microalgae model organism

Marco Di Marcello<sup>1</sup>, Iglesias-Gonzalez J<sup>1</sup>, Meszaros A<sup>1,2</sup>, Haider M<sup>3</sup>, Gnaiger E<sup>1,2</sup>, Huete-Ortega M<sup>1</sup>

<sup>1</sup>Oroboros Instruments, Innsbruck, Austria; <sup>2</sup>D Swarovski Research Lab, Dept Visceral, Transplant Thoracic Surgery, Med Univ Innsbruck, Austria; <sup>3</sup>SHTech, Innsbruck, Austria

✉ [marco.dimarcello@orooboros.at](mailto:marco.dimarcello@orooboros.at)

[http://www.mitofit.org/index.php/Di\\_Marcello\\_2019\\_MitoFit\\_Preprint\\_Arch](http://www.mitofit.org/index.php/Di_Marcello_2019_MitoFit_Preprint_Arch)



© 2019 Di Marcello *et al.* This is an Open Access extended abstract (not peer-reviewed) distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original authors and source are credited. © remains with the authors, who have granted MitoFit an Open Access preprint license in perpetuity.

**Editor MitoFit Preprint Archives:** Gnaiger E

### Introduction

Bioenergetics is the study of how living organisms acquire and transform energy to perform biological work. Energetic coupling between chloroplasts and mitochondria has been described in algae, demonstrating that a good functionality and interaction between both organelles is necessary to maintain metabolic integrity. High-resolution respirometry (HRR) is widely used to assess mitochondrial respiration and other bioenergetics parameters in the biomedical field of mitochondrial research and its clinical applications [1]. In our interdisciplinary study, we adapted the multimodal approach of the Oroboros O2k high-resolution respirometer to investigate algal bioenergetics for biotechnological purposes [2].

In contrast to mammalian cells, algal mitochondria possess alternative oxidases (AOX), which bypass electron transfer from the Q-junction through Complexes CIII and CIV [3]. Therefore, in algae we can distinguish between respiration through the Q-AOX and Q-CIV branches.


### Material and methods

The microalgal model organism *Chlamydomonas reinhardtii* wild-type strain *wt12* was grown at RT in Tris-Acetate-Phosphate (TAP) medium in a 16:8 h light:dark cycle. Oxygen flux,  $J_{O_2}$ , was monitored in *wt12* living cells in the exponential growth phase at 25 °C in Oroboros O2k high-resolution respirometers excluding any light in the chambers. Substrate-uncoupler-inhibitor titration (SUIT) protocols were specifically developed to characterise activities of the Q-AOX and Q-CIV branch (SUIT-022 [4] and SUIT-023 O<sub>2</sub> [5], respectively). To quantify the contribution of the Q-AOX branch to algal dark respiration, we studied the oxygen kinetics of (1) ROUTINE-respiration in TAP medium, (2) Q-AOX dependent respiration after inhibition of CIV with 1 mM potassium cyanide (KCN), and (3) Q-CIV dependent respiration after inhibition of AOX with 1 mM salicylhydroxamic acid (SHAM). Oxygen kinetics was obtained from aerobic-anaerobic


transitions with high time resolution at a data sampling interval of 0.2 s.  $p_{50}$  is the  $O_2$  partial pressure,  $p_{O_2}$ , at 50% of maximal respiration,  $J_{max}$  [6]. The  $p_{50}$  was calculated from hyperbolic fits using the Oroboros O2Kinetics software for automatic  $O_2$  calibration, correction for zero  $O_2$  signal drift, instrumental background  $O_2$  flux and exponential time constant of the polarographic oxygen sensor [7]. A single shifted hyperbolic fit was used to fit  $J_{O_2}$  as a function of  $p_{O_2}$  in each aerobic-anaerobic transition.

## Results and conclusions

$p_{50}$  ranged from 0.06 to 0.08 kPa for ROUTINE-respiration with an excellent fit by a first-order hyperbolic function. This oxygen affinity is comparable to that in small mammalian cells [8]. Upon inhibition of CIV with KCN,  $J_{O_2}$  was significantly impaired (Fig. 1A) and  $p_{50}$  increased three-fold up to 0.35 kPa (Fig. 2). No decline of  $J_{O_2}$  and  $p_{50}$  was observed relative to ROUTINE-respiration after inhibition of AOX with SHAM (Fig. 1B). In all cases, excellent fits of respiration as a function of oxygen pressure were obtained by a first-order hyperbolic function.

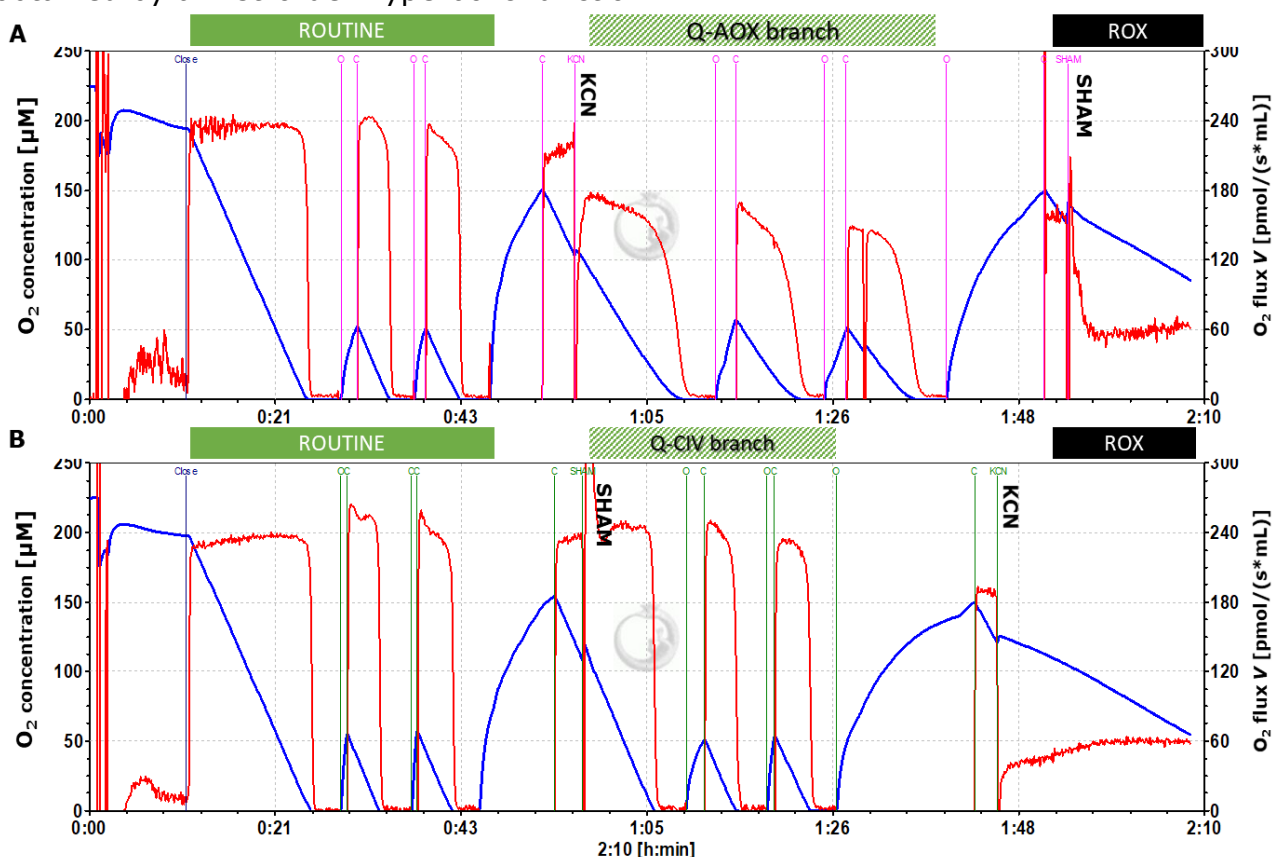

**Figure 1. High-resolution respirometry for the study of dark respiration and  $O_2$  kinetics with *C. reinhardtii* wt12.** Representative O2k traces showing  $O_2$  concentration and  $O_2$  flux per chamber volume with repeated aerobic-anoxic transitions ( $O_2$  kinetics) and re-oxygenations. **A:** Protocol SUIT-022: AOX-ce CN+SHAM. **B:** Protocol SUIT-023: AOX-ce SHAM+CN. Note the high technical reproducibility of ROUTINE-respiration in both protocols, and the identical and relatively high residual oxygen consumption,  $Rox$ , after titration of both inhibitors in both protocols.

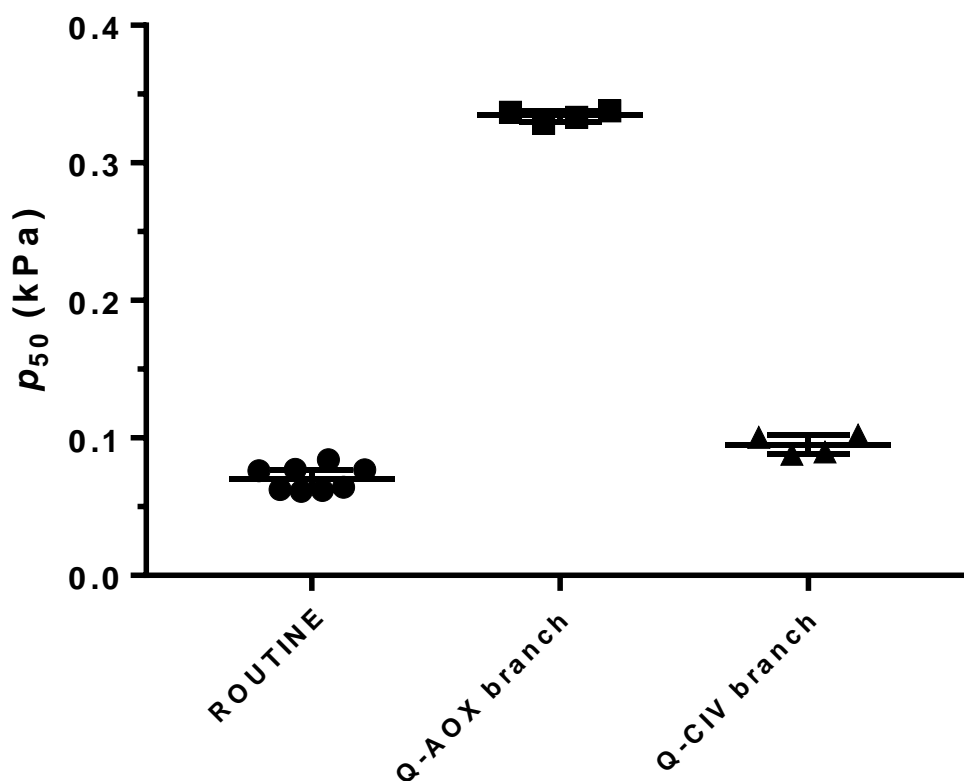

**Figure 2.  $p_{50}$  in living cells of *C. reinhardtii* in the ROUTINE-state of respiration, and metabolic pathways restricted to the Q-AOX or Q-CIV branch.**  $O_2$  kinetic experiments were run in presence of the cytochrome *c* oxidase inhibitor potassium cyanide (AOX group) or the alternative oxidase inhibitor salicylhydroxamic acid (CIV group). The data represents  $n=8$  technical replicates,  $N=2$ , median  $\pm$  IR.

If the potential contribution of the Q-AOX branch in the ROUTINE-state would be compensated for by increased Q-CIV flux after addition of SHAM, then the mixed Q-AOX and Q-CIV fluxes would give rise to biphasic hyperbolic oxygen kinetics, with a contribution of the high-affinity Q-CIV branch and the low-affinity Q-AOX branch. Taken together, our results provide evidence against a contribution of AOX to ROUTINE-dark respiration in *wt12* cells under the presently applied culture conditions. Oxygen kinetics provides a sensitive and fast method for detection of Q-AOX and Q-CIV contributions to dark respiration in living cells. This kinetic approach is based on the difference of  $O_2$  affinities of the two pathway branches, which will extend our understanding of the bioenergetics and physiology of all types of cells harbouring AOX and CIV.



Supported by project NextGen-O2k, which has received funding from the European Union's Horizon 2020 research and innovation programme under grant agreement No 859770

## References

1. Doerrier C, Garcia-Souza LF, Krumschnabel G, Wohlfarter Y, Mészáros AT, Gnaiger E (2018) High-Resolution Fluorescence Respirometry and OXPHOS protocols for human cells, permeabilized fibers from small biopsies of muscle, and isolated mitochondria. *Methods Mol Biol* 1782:31-70. - [www.bioblast.at/index.php/Doerrier\\_2018\\_Methods\\_Mol\\_Biol](http://www.bioblast.at/index.php/Doerrier_2018_Methods_Mol_Biol)

2. Pulz O & Gross W (2004) Valuable products from biotechnology of microalgae. *App. Microbiol. Biotech.* 65(6): 635-648.
3. Young L, Shiba T, Harada S, Kita K, Albury MS, Moore AL (2013) The alternative oxidases: simple oxidoreductase proteins with complex functions. *Biochem Soc Trans* 41:1305-11.
4. [www.bioblast.at/index.php/SUIT-022\\_O2\\_ce\\_D051](http://www.bioblast.at/index.php/SUIT-022_O2_ce_D051)
5. [www.bioblast.at/index.php/SUIT-023\\_O2\\_ce\\_D053](http://www.bioblast.at/index.php/SUIT-023_O2_ce_D053)
6. Gnaiger E (2001) Bioenergetics at low oxygen: dependence of respiration and phosphorylation on oxygen and adenosine diphosphate supply. *Respir Physiol* 128:277-97. - [www.bioblast.at/index.php/Gnaiger\\_2001\\_Respir\\_Physiol](http://www.bioblast.at/index.php/Gnaiger_2001_Respir_Physiol)
7. Meszaros AT, Haider M, Di Marcello M, Gnaiger E (2018) High-resolution mitochondrial oxygen kinetics as diagnostic tool in Complex IV impairments. Abstract Mitochondrial Medicine 2018 Hinxton UK. - [www.bioblast.at/index.php/Meszaros\\_2018\\_Mt\\_Med\\_Hinxton](http://www.bioblast.at/index.php/Meszaros_2018_Mt_Med_Hinxton)
8. Scandurra FM, Gnaiger E (2010) Cell respiration under hypoxia: facts and artefacts in mitochondrial oxygen kinetics. *Adv Exp Med Biol* 662:7-25. - [www.bioblast.at/index.php/Scandurra\\_2010\\_Adv\\_Exp\\_Med\\_Biol](http://www.bioblast.at/index.php/Scandurra_2010_Adv_Exp_Med_Biol)